

The normal electroretinogram in adult healthy Shih Tzu dogs using the HMsERG

June-sub Lee¹, Kyung-hee Kim¹, Ha-young Jang¹, Bora Lee¹, Joon Young Kim², Soon-wuk Jeong^{1,*}

¹Department of Veterinary Surgery, College of Veterinary Medicine, Konkuk University, Seoul 143-701, Korea

²Department of Veterinary Medicine, Madingley Road, Cambridge CB3 0ES, United Kingdom

Electroretinography (ERG) is a reliable diagnostic tool for the diagnosis of retinal disease. It measures electric potentials occurring in the retina in response to light stimulation. In this study, we examined the normal electroretinogram using the Handheld Multispecies ERG (HMsERG) in Shih Tzu dogs. ERG recordings were performed in twelve eyes of six healthy Shih Tzu dogs. Dogs were anesthetized with a combination of medetomidine and ketamine. Proparacaine eye drops were also applied as a topical anesthetic. Tropicamide eye drops were applied for mydriasis. After 20 min of dark adaptation, we recorded the amplitudes and implicit times of the b-waves of the rod, standard rod and cone (Std R&C), high-intensity rod and cone (Hi-int R&C), and cone systems, and responses of the cones and inner retina by flicker light stimulation (cone flicker). Results showed that mean the amplitudes of a-waves of Std R&C, Hi-int R&C, and the cone responses were 141.25 μ V, 173.00 μ V, and 12.92 μ V, respectively. The b-waves of the rod responses ranged from 141.58 to 155.25 μ V; the Std R&C was 314.75 μ V, the Hi-int R&C was 329.42 μ V, the cones were 37.75 μ V, and the flicker responses were 64.08 μ V. The b/a ratios for the Std R&C, Hi-int R&C, and the cone response were 2.29, 1.94, and 3.71, respectively. Mean implicit time of the a-wave of the Std R&C was 15.12 ms, of Hi-int R&C was 13.42 ms, and of the cone response was 7.22 ms. The b-wave of the rod responses ranged from 68.12 to 72.68 ms, of Std R&C were 37.28 ms, of Hi-int R&C were 41.90, of the cone responses were 38.12 ms, and of the cone flicker responses were 22.80 ms. We believe that these parameters can be used as reference “normal” ERGs ranges for Shih Tzu dogs using the HMsERG under medetomidine and ketamine anesthesia.

Keywords: electroretinogram, HMsERG, Shih Tzu dog

Introduction

Electroretinography (ERG) is a reliable diagnostic tool for the evaluation of retinal function. It is commonly used to diagnose retinal disorders before cataract extraction or in cases of unexplained visual loss without ophthalmoscopic abnormalities [11].

The definition of a normal range of ERG values is difficult because the results of ERG can vary due to intrinsic factors and extrinsic factors. Of the former, the eye's state of light adaptation affects the results of the ERG recording the most, followed by species, age, transparency of the ocular media, retinal integrity, retinal circulatory disturbances, ocular opacity and pupil dilation [2]. Extrinsic factors are time and intensity for light stimulation, location and type of electrodes, kind of recording equipment used, anesthetic protocol, experimental conditions, and environment factors [2]. All of these factors may differ between various laboratories. However, yielding normal ranges for the ERG will certainly help to diagnose retinal disease processes. Normal ranges should be obtained for each specific clinic or laboratory, breed, age range, and every ERG equipment used.

In this study, the amplitude and implicit time ranges of the a- and b- wave, and b/a ratio of normal eyes in healthy Shih Tzu dogs were measured using the Handheld Multispecies ERG (HMsERG) for determining the standard ERG parameters in our laboratory.

Materials and Methods

Animals

Twelve eyes of six healthy Shih Tzu dogs were used. All were normal based on the physical, hematological, serological and ocular examination results. Their average age and weight were 2.28 ± 0.48 (mean \pm SE) and 3.59 ± 6.45 kg, respectively. They were individually confined and fed commercial dry food and water *ad libitum*.

*Corresponding author

Tel: +82-2-450-3670; Fax: +82-2-456-4655

E-mail: swjeong@konkuk.ac.kr

ERG equipment

The portable mini-Ganzfeld ERG unit (HM_sERG; RetVet, USA) was used in this study. The frequency band was 0.3 to 300 Hz. The primary components of the ERG unit were the stimulation and the recording systems. The unit contained a compact flash card (CFC) with software to drive the ERG unit. Obtained ERG recordings were stored on the CFC. The reference and ground electrodes were needle electrodes (Model E2-straight needle; Astro-Med, USA) and an ERG Jet contact lens electrode (Nicolet Biomedicals, USA) was used as the active electrode. The obtained responses were transferred to a computer to visualize the ERG readings and to print the results.

ERG procedure

ERG examinations were performed under dim red light in a dark room. The pupils were dilated with 0.5% tropicamide (eye drops, two times, 15 min intervals; SAMIL, Korea). After producing mydriasis, animals were dark-adapted for 20 min [4,10,12,15]. ERG examination was then performed under general anesthesia using a combination of 40 µg/kg medetomidine hydrochloride (Orion Pharma, Finland) and 4 mg/kg ketamine hydrochloride (Yuhan, Korea) intravenously [6,8,12]. The ground electrode was placed over the external occipital protuberance and the reference electrode was placed approximately midpoint between lateral canthus and ipsilateral ear pinna [8]. The active electrode was positioned on the cornea after topical anesthetic eye drops (0.5% proparacaine hydrochloride ophthalmic solution; Alcon, Belgium) and artificial tears (0.03% hydroxypropylmethyl cellulose; Hanlim, Korea) were applied [2,3,14]. Before ERGs were recorded, impedance and baseline tests were performed; the latter for evaluating the noise level in the environment. ERGs were recorded automatically using the dog diagnostic protocol pre-programmed for this equipment which came from the manufacturer that way. It consists of a total of 17 sessions and the entire protocol takes 34.57 min per eye. All functions, including intensity of light and time for each session are described in Table 1 [5].

Evaluation of ERG

The a-wave amplitude was measured from the baseline to the a-wave trough, and the b-wave amplitude was measured from the a-wave trough to the b-wave peak [9]. Then b/a ratios of amplitude were calculated. The a- and b-wave implicit times were measured from the stimulus onset to the a-wave trough and b-wave peak, respectively [9]. For estimation of normal ranges of each measurement values, the mean and standard error of a- and b-wave amplitudes, and implicit times, as well as b/a ratios were calculated. Because there is no a-wave in the dark adapted eye at 10 mcd.s/m², the a-wave at 10 mcd.s/m² were not used.

Table 1. Functions, intensity of light and time for each electroretinography session

Session	Function	Intensity of light (mcd.s/m ²)	Time (second)
1st	Dark adaptation cycle		240
2nd	Rod response	10	20
3rd	Dark adaptation cycle		240
4th	Rod response	10	20
5th	Dark adaptation cycle		240
6th	Rod response	10	20
7th	Dark adaptation cycle		240
8th	Rod response	10	20
9th	Dark adaptation cycle		240
10th	Rod response	10	20
11th	Delay		2
12th	Std R&C response	3,000	40
13th	Delay		30
14th	Hi-int R&C response	10,000	80
15th-1	Light adaptation cycle	3,000	600
15th-2	Cone response	3,000	16
16th	Delay		2
17th	Cone flicker response	3,000	4.1

Std R&C: standard rod and cone, Hi-int R&C: high-intensity rod and cone.

Results

The results of ERG recordings using the dog diagnostic protocol of the HM_sERG are represented as waveforms (Fig. 1). Both a- and b- waves were represented in the standard rod and cone (Std R&C), high-intensity rod and cone (Hi-int R&C), and cone responses, but only b-waves were obtained for the rod and cone flicker responses (Fig. 1).

a- and b-wave amplitude and b/a ratio

The a-wave amplitude is shown in Fig. 2 and Table 2. The Hi-int R&C response was high at 173 ± 8.80 µV (mean ± SE); the Std R&C and the following cone responses were 141.25 ± 9.04 µV and 12.92 ± 3.97 µV, respectively.

The b-wave amplitude is shown in Fig. 3 and Table 2. Rods 1 to 5 had amplitudes of 150.08 ± 14.36 µV, 145.50 ± 14.69 µV, 141.92 ± 16.06 µV, 141.58 ± 18.55 µV and 155.25 ± 14.60 µV, respectively. It was increased in the Std R&C response with 314.75 ± 17.97 µV and the highest in the Hi-int R&C response with 329.42 ± 16.47 µV. The cone and cone flicker responses had lower amplitudes than the rod responses at 37.75 ± 4.90 µV and 64.08 ± 5.26 µV, respectively.

The b/a ratios for Std R&C, Hi-int R&C and cone responses were 2.29 ± 0.15 µV, 1.94 ± 0.12 µV, and 3.71 ± 0.47 µV, respectively (Fig. 4 and Table 2).

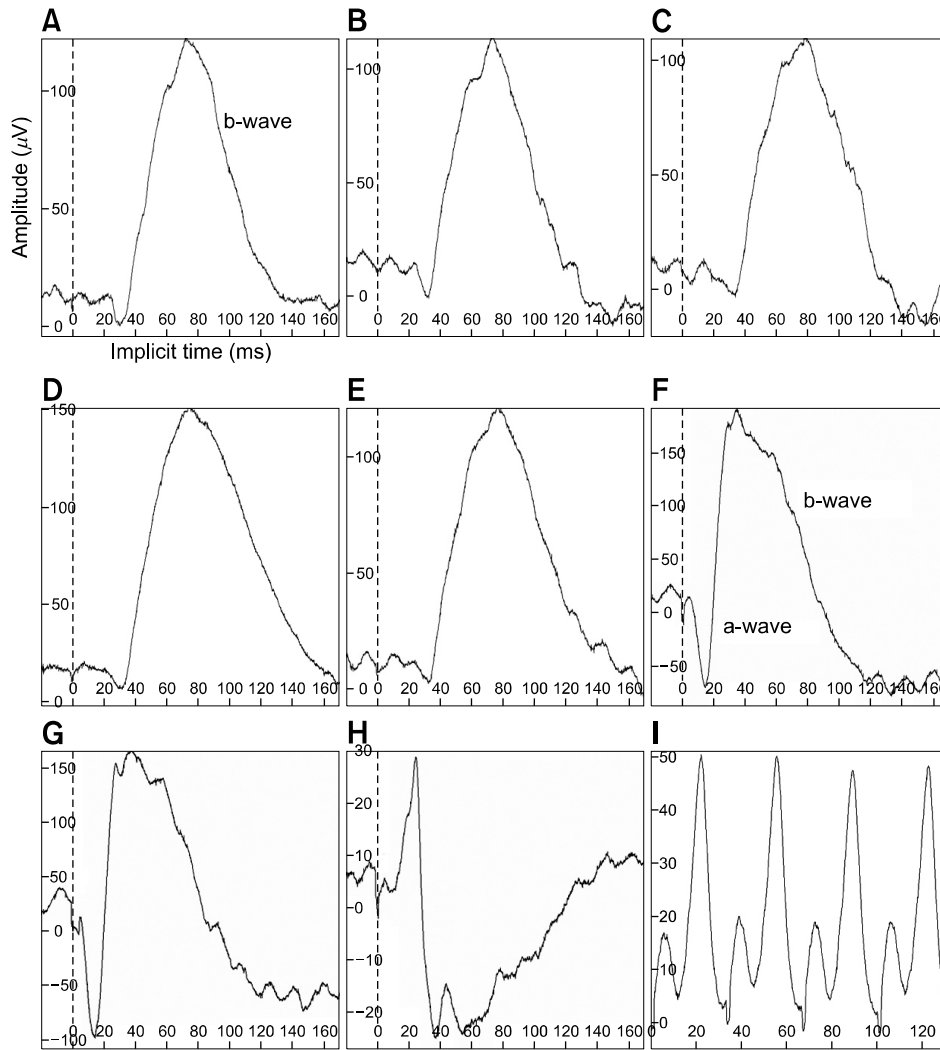


Fig. 1. Representative waveforms of the ERG results using the dog diagnostic protocol of the Handheld Multispecies electroretinography (HMsERG) in adult Shih Tzu dogs showing normal retinal function. (A-E: rod response, F: standard rod and cone (Std R&C) responses, G: high-intensity (Hi-int) R&C responses, H: cone response, I: cone flicker response).

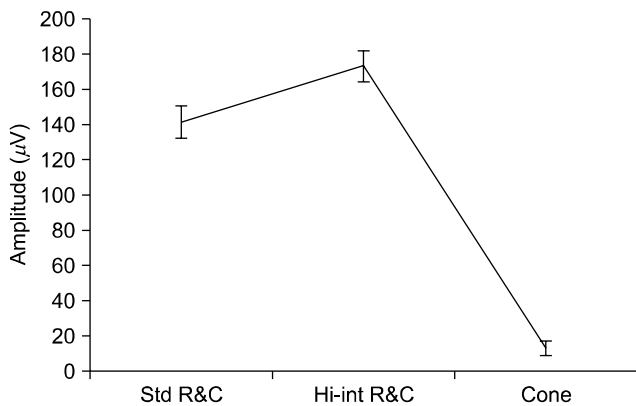


Fig. 2. a-wave amplitudes (mean ± SE) using the diagnostic protocol with the HMsERG equipment in adult Shih Tzu dogs showing normal vision. Std R&C: standard rod and cone, Hi-int R&C: high-intensity rod and cone.

a- and b-wave implicit time

The a-wave implicit time is shown in Fig. 5 and Table 3. It was the fastest for the Std R&C response at 15.12 ± 0.17 ms (mean ± SE), and were 13.42 ± 0.39 ms and 7.22 ± 3.23 ms for the Hi-int R&C and cone responses, respectively.

The b-wave implicit time is shown in Fig. 6 and Table 3. The implicit time of rods 1-5 were the highest at 72.58 ± 1.85 ms, 70.86 ± 2.17 ms, 70.65 ± 2.32 ms, 68.12 ± 5.38 ms and 72.68 ± 1.83 ms, respectively. Std R&C and Hi-int R&C responses showed similar values at 37.28 ± 2.67 ms and 41.90 ± 2.03 ms, respectively, and then it decreased for the cone response to 38.12 ± 12.13 ms.

Discussion

Shih Tzus were selected for this study because it is the

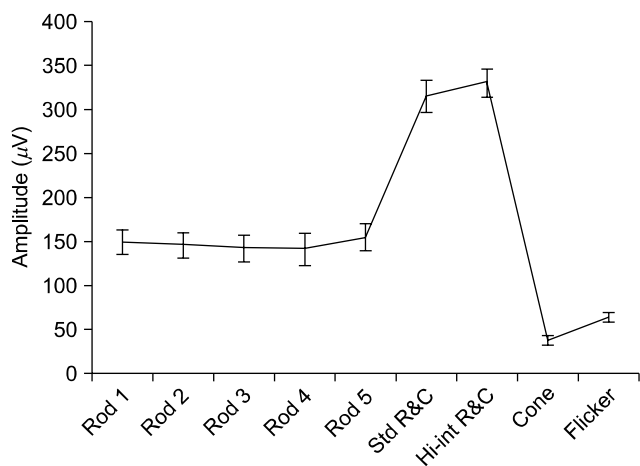


Fig. 3. b-wave amplitudes (mean ± SE) using the dog diagnostic protocol with the HM_sERG equipment in adult Shih Tzu dogs showing normal vision. Std R&C: standard rod and cone, Hi-int R&C: high-intensity rod and cone.

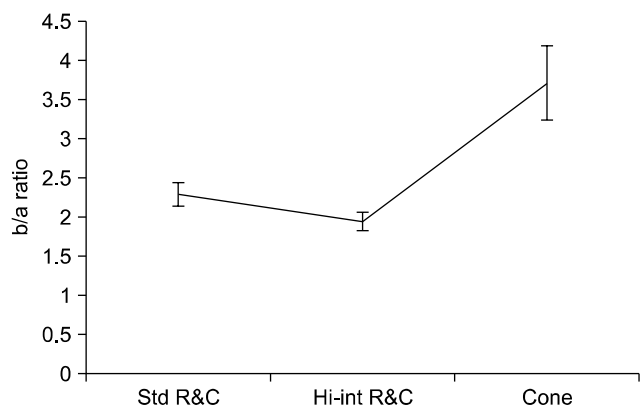


Fig. 4. b/a ratio of amplitudes (mean ± SE) using dog diagnostic protocol of HM_sERG equipment in adult Shih Tzu dogs showing normal vision. Std R&C: standard rod and cone, Hi-int R&C: high-intensity rod and cone.

most preferred and the most common dog breed in Korea. The breed has a predisposition for eye diseases which requires retinal evaluation for diagnosis, such as retinal detachment, progressive retinal atrophy and cataracts [13].

Various anesthetic protocols can affect the results of ERG examinations. According to previous studies, halothane and sevoflurane strongly depress the scotopic threshold response in Beagles, while moderately depressing the b-wave and increasing the oscillatory potential amplitudes, precluding their use for ERG examinations [14]. On the other hand, the combination of ketamine and xylazine has been commonly used because the ketamine and xylazine mix achieves sufficient immobilization for disturbance-free ERG recordings [3]. In regards to ketamine, Kommonen and Raitta [3] reported the most adequate time for the ERG recording was 10 to 20 min after ketamine administration.

Table 2. a- and b- wave amplitude and b/a ratio

Response	a-wave amplitude (μV)	b- wave amplitude (μV)	b/a ratio
Rod 1	None	150.08 ± 14.36 (68 to 248)	None
Rod 2	None	145.50 ± 14.69 (61 to 240)	None
Rod 3	None	141.92 ± 16.06 (48 to 234)	None
Rod 4	None	141.58 ± 18.55 (7 to 229)	None
Rod 5	None	155.25 ± 14.60 (82 to 231)	None
Std R&C	141.25 ± 9.04 (89 to 211)	314.75 ± 17.97 (245 to 474)	2.29 ± 0.15 (1.61 to 3.33)
Hi-int R&C	173 ± 8.8 (132 to 240)	329.42 ± 16.47 (255 to 460)	1.94 ± 0.12 (1.44 to 2.86)
Cone	12.92 ± 3.97 (6 to 56)	37.75 ± 4.90 (23 to 72)	3.71 ± 0.47 (1.29 to 6.8)
Cone flicker	None	64.08 ± 5.26 (42 to 103)	None

Std R&C: standard rod and cone, Hi-int R&C: high-intensity rod and cone. Mean ± SE (min. to max.).

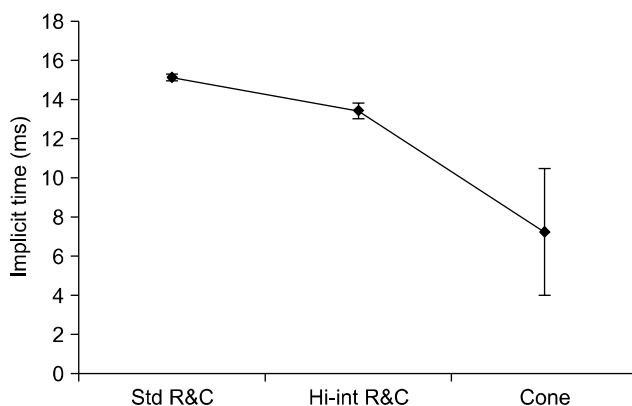


Fig. 5. a-wave implicit time (mean ± SE) using dog diagnostic protocol of HM_sERG equipment in adult Shih Tzu dogs showing normal vision. Std R&C: standard rod and cone, Hi-int R&C: high-intensity rod and cone.

Furthermore, it had no significant effect on the ERG results. This protocol also showed that there was no globe rotation for Poodle and Beagle dogs. Therefore it did not require stay-sutures for restriction of eyeball rotation. However, this effect cannot be achieved in Labrador Retrievers [3]. Isoflurane appears not to affect the results of ERG recordings. But eyeball rotation may occur, so stay-sutures are recommended [4,5].

In this study, ketamine and medetomidine was used intravenously. We used medetomidine instead of xylazine because of the vomiting risk [5]. With this combination, appropriate anesthesia was achieved for ERG examination, which took almost 40 min. Moreover, eyeball rotation was not observed in any of the dogs. Because atipamezol

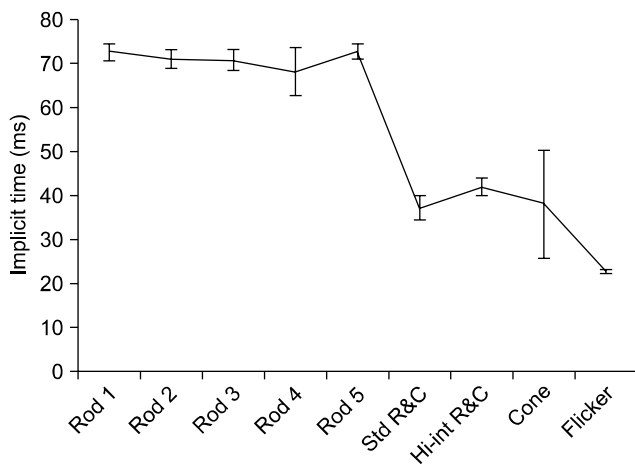


Fig. 6. b-wave implicit time (mean \pm SE) using dog diagnostic protocol of HM sERG equipment in adult Shih Tzu dogs showing normal vision. Std R&C: standard rod and cone, Hi-int R&C: high-intensity rod and cone.

hydrochloride could be applied as a counter if needed, we believe that this protocol is safer than previous anesthetic protocols. Therefore, the combination of ketamine and medetomidine can be considered as a good anesthetic protocol to use with the HM sERG in the Shih Tzu dog.

The HM sERG consists of the unit body, a mini-ganzfeld and three electrodes; the active electrode (the jet lens electrode) and 2 needle electrodes. Very advantageous in regards to the HM sERG is the mini-ganzfeld. Up to now, the conventional Ganzfeld dome has been widely used for full-field ERG examination in veterinary ophthalmology [5,7,9]. However, some disadvantages of the conventional Ganzfeld dome are the large and expensive system, the difficulty of using bilateral stimulation recordings, and the fixation of the animal during examination [5]. Also the location of the dog's head inside a sphere of 60 cm diameter may be difficult to achieve [5]. But the mini-ganzfeld overcomes these issues.

The HM sERG has four built in protocols based on the International Society for Clinical Electrophysiology of Vision (ISCEV). Among them, the dog diagnostic protocol was used in this study because this provides more accurate values of rod and cone responses in dogs than the case for short protocols [1].

Generally, ERG examinations are interpreted with amplitude and implicit times of a- and b-waves. In addition, one of the reasons that the evaluation of b/a ratios are considered to be an important parameter being that it is an indicator of disorders of the retina in cases with dense opacity observed in the anterior segment and vitreous body [5].

The ISCEV established standard flash b/a ratio = 2 at intensity of 3 cds./m² in humans [7]. Maehara *et al.* [5] reported similar b/a ratio at the same intensity of light in beagles. In this study, b/a ratio was 2.29 \pm 0.15 at the intensity

Table 3. a - and b - wave implicit time

Response	a-wave implicit time (ms)	b-wave implicit time (ms)
Rod 1	None	72.58 \pm 1.85 (60.3 to 84.9)
Rod 2	None	70.86 \pm 2.17 (54 to 82.5)
Rod 3	None	70.65 \pm 2.32 (59 to 83.2)
Rod 4	None	68.12 \pm 5.38 (10.6 to 83.5)
Rod 5	None	72.68 \pm 1.83 (63.3 to 82.6)
Std R&C	15.12 \pm 0.17 (14 to 16)	37.28 \pm 2.67 (30.9 to 59)
Hi-int R&C	13.42 \pm 0.39 (11.9 to 15.9)	41.90 \pm 2.03 (30.4 to 53.4)
Cone	7.22 \pm 3.23 (0.6 to 39.2)	38.12 \pm 12.13 (24.8 to 171.5)
Cone flicker	None	22.80 \pm 0.63 (16.2 to 24.5)

Std R&C: standard rod and cone, Hi-int R&C: high-intensity rod and cone. Mean \pm SE (min. to max.).

3 cds./m². This result is similar to the above studies.

The time of dark adaptation is important since it can affect the ERG results as the values may be affected by the degree of light exposure before the ERG [11]. Because the dog diagnostic protocol has a dark adaptation cycle (a low intensity stimulation for every 4 min) of 20 min, any further dark adaptation is not needed before performing ERG examinations. However, when we applied 20 min of dark adaptation before performing an ERG examination, then we obtained the same b-wave amplitudes as the dark adaptation cycle. This means we were not checking the dark adaptation cycle properly. But these results can be used as the normal values of specific responses such as the rod system b-wave, the mixed response rod and cone a- and b- waves (for scotopic recordings) and the cone a- and b-waves and the flicker b-waves (for photopic recordings).

In conclusion, we believe that the parameters obtained in this study can be used as normal ERG reference ranges for Shih Tzu dogs.

References

1. Ford M, Bragadóttir R, Rakoczy PE, Narfström K. Gene transfer in the RPE65 null mutation dog: relationship between construct volume, visual behavior and electroretinographic (ERG) results. *Doc Ophthalmol* 2003, **107**, 79-86.
2. Gum GG. Electrophysiology in veterinary ophthalmology. *Vet Clin North Am Small Anim Pract* 1980, **10**, 437-454.
3. Kommonen B, Raitta C. Electroretinography in Labrador Retrievers given ketamine-xylazine anesthesia. *Am J Vet Res* 1987, **48**, 1325-1331.
4. Krill AE. The electroretinogram and electro-oculogram: clinical applications. *Invest Ophthalmol* 1970, **9**, 600-617.
5. Maehara S, Itoh N, Itoh Y, Wakaiki S, Tsuzuki K, Seno

- T, Kushiro T, Yamashita K, Izumisawa Y, Kotani T.** Electroretinography using Contact Lens Electrode with Built-in Light Source in dogs. *J Vet Med Sci* 2005, **67**, 509-514.
6. **Marmor MF, Holder GE, Seeliger MW, Yamamoto S.** Standard for clinical electroretinography (2004 update). *Doc Ophthalmol* 2004, **108**, 107-114.
 7. **Marmor MF, Zrenner E.** Standard for clinical electroretinography (1999 update). International Society for Clinical Electrophysiology of Vision. *Doc Ophthalmol* 1999, **97**, 143-156.
 8. **Mentzer AE, Eifler DM, Montiani-Ferreira F, Tuntivanich N, Forcier JQ, Petersen-Jones SM.** Influence of recording electrode type and reference electrode position on the canine electroretinogram. *Doc Ophthalmol* 2005, **111**, 95-106.
 9. **Narfström K, Ekesten B, Rosolen SG, Spiess BM, Percicot CL, Ofri R.** Guidelines for clinical electroretinography in the dog. *Doc Ophthalmol* 2002, **105**, 83-92.
 10. **Sims MH, Brooks DE.** Changes in oscillatory potentials in the canine electroretinogram during dark adaptation. *Am J Vet Res* 1990, **51**, 1580-1586.
 11. **Slatter DH.** *Fundamentals of Veterinary Ophthalmology*. 3rd ed. pp. 419-456, Saunders, Philadelphia, 2001.
 12. **Tuntivanich N, Mentzer AL, Eifler DM, Montiani-Ferreira F, Forcier JQ, Johnson CA, Petersen-Jones SM.** Assessment of the dark-adaptation time required for recovery of electroretinographic responses in dogs after fundus photography and indirect ophthalmoscopy. *Am J Vet Res* 2005, **66**, 1798-1804.
 13. **Whitley RD, Vygantas KR.** Presumed inherited ocular diseases. In: Martin CL (ed.). *Ophthalmic Disease in Veterinary Medicine*. 1st ed. pp. 471-491, Manson, London, 2005.
 14. **Yanase J, Ogawa H.** Effects of halothane and sevoflurane on the electroretinogram of dogs. *Am J Vet Res* 1997, **58**, 904-909.
 15. **Yanase J, Ogawa H, Ohtsuka H.** Scotopic threshold response of the electroretinogram of dogs. *Am J Vet Res* 1996, **57**, 361-366.